# happylegs


**Instrucciones**  
**User Manual**  
**Bedienungsanleitung**  
**Mode d'emploi**

[www.happylegs.es](http://www.happylegs.es)

# happylegs®

la máquina de andar sentado

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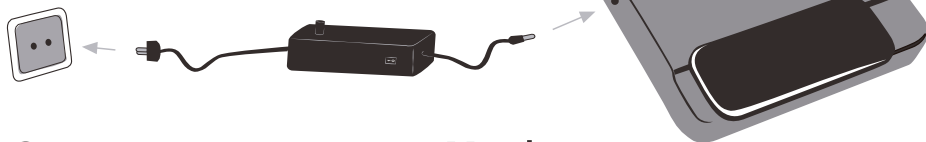
the seated walking machine

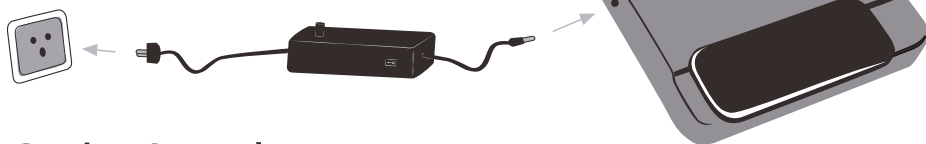
## Getting Started.

### The box contains 1 Happylegs and 1 Power Laptop Transformer:

The powered exerciser Happylegs has an energy-efficient motor of low power consumption. The equipment connects to the mains through three-speed Laptop transformer. To use Happylegs simply connect the transformer to the connector at the front of the Happylegs machine and then plug the transformer socket into mains electricity. Turn the machine on (1) and off (0) from the transformer switch and turn the dial to select the required speed.

If you want to lean the Happylegs unscrewing the two height-adjustable legs in the lower side. Sit back, put your feet on the pedals, relax and let Happylegs move your feet and legs with its smooth motion. Do not force the motion, because Happylegs automatically moves your legs without any effort.
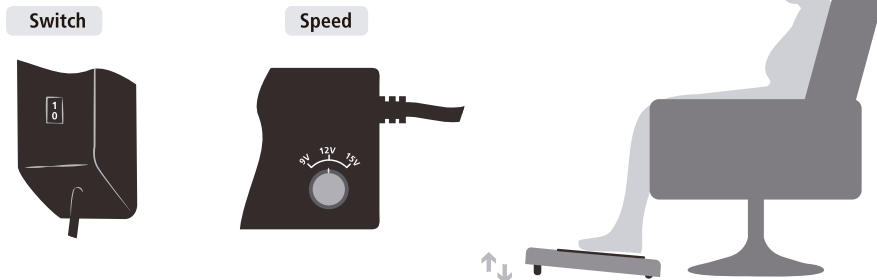
## Select the required operational speed.

The speed can be altered by selected differing voltages with the speed selector:

- **9 V Slowest:** This speed is the more gentle an ideal for people who may be recovering for injury or have stiff joints and also for those who suffer fatigue.
- **12 V Normal Speed:** Recommended for prolonged use at home or at work. This should produce a pleasant and relaxed feeling.
- **15 V Fastest Speed:** The use of this speed is recommended for periods of less than two hours continuous to elder people and more prolonged to middle-aged.

## **Cleaning and Maintenance.**

In order to use Happylegs correctly and for safety reasons, following the instructions outlined should be followed. Never stand on top of the machine, as full weight could damage the internal operational mechanism. The machine should not be submitted to pressure for which is not designed. Happylegs does not require regular maintenance and has been developed to be used continuously for long periods of time. In order to clean it, a dry cloth should be used when the equipment is unplugged and switched off.

## **Transportation and Storage.**

This product is electrically operated avoid spilling any liquid on the Happylegs or staying in high humidity environments. Also avoid temperatures greater than 60° C, which could damage the equipment. Happylegs should be transported within its packaging. Avoid bumping or dropping and do not store weight greater than 40 kg on it. If broken or damaged contact your official supplier. Not to be disassembled or repaired by unauthorized personnel.

## **Technical Specifications.**

Happylegs is powered by mains electricity via a transformer designed exclusively for Autoergon. Do not use parts or accessories not supplied by the official Happylegs manufacturer or supplier.

- Connected to the mains electricity.
- Input Voltage: 110-240 V// 50-60Hz 1,5A 
- Output Voltage: 9V/ 12V/ 15V 3A 

B.  Class II     RoHS

## **Warranty Terms and CE Marking.**

Happylegs has been manufactured with the highest quality both in the materials used and quality control and design of its components. This has led to passing the most stringent quality tests needed to obtain the European Quality certification (CE Marking).


The Happylegs machine has a two-year warranty against all manufacturing defects. The warranty does not cover damage to the machine caused by dropping it or being knocked, nor does it cover any breakage to the power cable resulting from misuse or negligence. The three-speed Laptop transformer with connector is covered by a one-year warranty. Never pull the power cable to unplug the machine as this can damage both the plug and the socket. Should any component fail to work, please contact the official after-sales service department where you purchased the machine.

## **How to dispose of old and unwanted electrical and electronic appliances:**

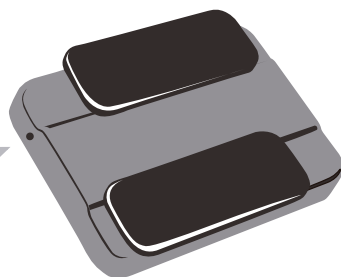1. If the symbol with a crossed container is displayed, this means it complies with Directive 2002/96/CE.
2. All electrical or electronic appliances should be disposed of properly at designated located places only.
3. The correct disposal and treatment of damaged goods contributes to avoiding potential risks to the public and the environment.
4. To obtain more information on how to properly dispose of your electrical and/or electronic devices, contact your local authorities or the retail place of purchase.

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